
Digital Bathroom Scale


Model HD-366

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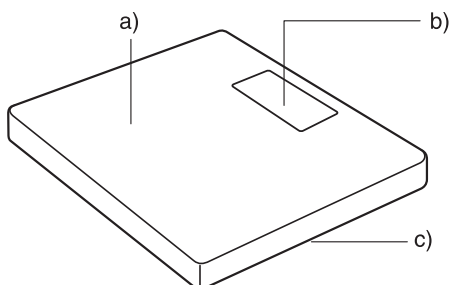2. (Fig.2)

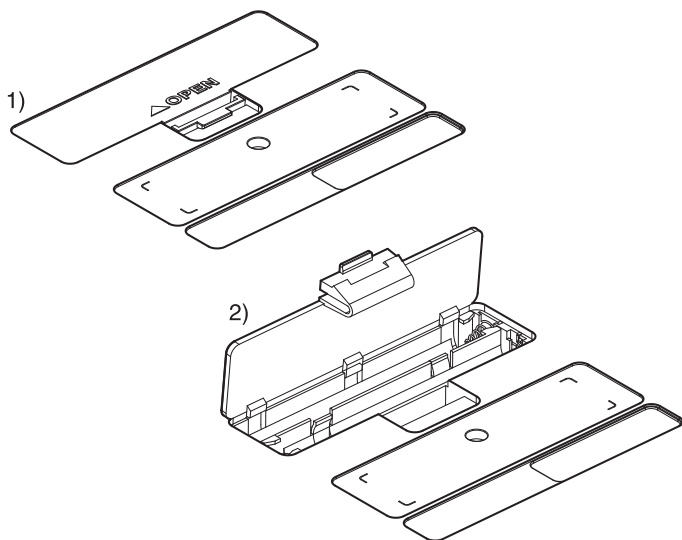**Eng**

This device features radio interference suppression in compliance with EC Regulation 2004/108/EC.

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2. Parts (Fig.1)

- a) Body
- b) Display
- c) Battery (Bottom)

3. Activating The Scale (Only When Installing New Batteries)

- 1) Insert batteries.
- 2) Place the scale on a firm even surface, making sure all four feet make contact at the same time.
- 3) With your foot tap the middle of the platform to activate the scale.
On the LCD display "8888" will appear – wait until it changes to "0.0" before stepping onto the platform.

4. Weighing

- 1) Select your preferred weight mode using the weight Mode Switch on the bottom of the scale.
- 2) Step onto the platform and stand still.
- 3) After few seconds your weight will flash on the display.
- 4) The scale will turn off automatically after 30 seconds.
 - * The scale should be used on a firm surface for accurate weighing, not recommended for use on carpeted surfaces.
 - * The LCD will display "OL", if the weight is over 200 kg / 440 lb / 31 st. 6 lb, and will turn off after 30 seconds.
 - * Whenever moving the scale to another location make sure that it is turned off.

5. Replacing the Batteries

When the low battery sign (Lo) appears on the LCD, remove the current batteries and replace them with a fresh set as indicated in Fig. 2.

After replacing batteries, repeat step #3 (Activating the scale).

6. Care Instructions

- 1) Be careful not to open the scale or tamper with its mechanism or circuitry.
These actions will void all warranty coverage.
- 2) Take care not to spill water or chemical liquids onto the scale.
Wipe the scale with a damp cloth periodically.
- 3) When not used for long periods of time, remove the batteries to prevent damage.
- 4) Do not use this scale while operating transmitters (such as a cell phone), as this may affect your results.
- 5) Do not step on the platform when wet.
- 6) Do not drop any objects onto the platform.

Federal Communications Commission and Canadian ICES Notice

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules and Canadian ICES-003.

These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio or television technician for help.

Modifications

The FCC requires the user to be notified that any changes or modifications made to this device that are not expressly approved by Tanita Corporation may void the user's authority to operate the equipment.

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