TANITA

Monitoring Your Health

Solar Scale

Balance solaire
Báscula (Pesa) Solar
Balança solar

Models 1618, 1631, 1632, HS-201, HS-301

Operating Instructions / Mode d'emploi / Instrucciones para la Operación / Instruções operacionais

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Federal Communications Commission Notice

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio or television technician for help.

Modifications

The FCC requires the user to be notified that any changes or modifications made to this device that are not expressly approved by Tanita Corporation may void the user's authority to operate the equipment.

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Features

- This scale is solar powered and never requires batteries.
- No complicated operation or set-up is required.

Before Use

Select your preferred weight mode (lb or kg) using the Weight Mode Switch.

Place your Tanita Solar Scale on a flat, hard surface with **direct exposure to natural or artificial lighting**. Tap the Reset Switch located on the lower right-hand edge of the measuring platform so that "0.0" appears on the LCD display. If the light source is too dim, no reading or a "Lo" message may appear, and you will need to move the scale to a brighter location.

Placing scale

The scale should be used on a hard even surface for accurate weighing. Not recommended for use on carpeted surfaces.

Instructions for Use


When using the scale for the first time, or after it has been moved, tap the Reset Switch with your toe so that "0.0" shows on the display (Fig. 1). Just step on the scale to measure your weight.

1. If the display shows other than "0.0," make sure the scale is exposed to a sufficient light source and tap the Reset Switch.
2. Step on the scale. The display shows body weight (Fig. 2).
3. Step off the scale. The display returns to "0.0" (Fig. 3).

Note: "OL" will appear on the display when the scale's weight capacity is exceeded. Avoid blocking light source when stepping onto scale.

Precautions

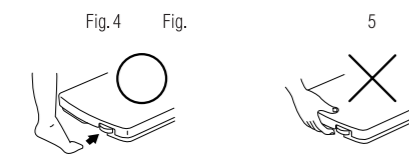To recalibrate the scale, gently tap the Reset Switch with your toe (see Fig. 4). Do not place your hand on the scale while recalibrating (see Fig. 5); if you do, "-0.5" or "-1.0" will be displayed. Using your toe, tap the Reset Switch again to display "0.0."

Care and Maintenance

- Place the scale in an area free from heating equipment, high humidity, or extreme temperature changes.
- Never submerge the scale in water or spill chemical liquids on it. To clean, use a dampened cloth and mild detergent.
- Avoid excessive impact or vibration to the unit.

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