TANITA®

BF-679W Body Fat Monitor

en Quick Reference Guide

- Fast, easy set up of your body fat monitor/scale.
- Refer to the Instruction Manual for more detailed information.

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General Programming

The unit will prompt you through the set up.
The display will highlight the function and the unit will beep after each step has been completed.

Use the Up/Down buttons to scroll through data.

Use the Set button to store data.

Set and Store Data


Turn on Power

Use the Set button to activate the unit.


Select Personal Data Number

Press the Up/Down buttons to select a personal key.
Press the Set button.

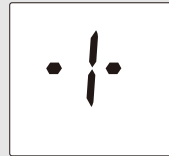
Set Age

Use the Up/Down buttons to select age.
Press the Set button.

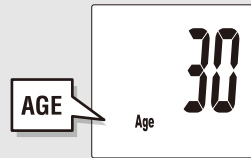
Select Female or Male

Use the Up/Down buttons to select your gender.
Press the Set button.

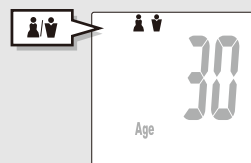Specify Height

Use the Up/Down buttons to select height.
Press the Set button.

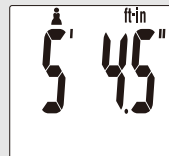The unit will beep twice and the display will flash all data three times to confirm the programming. The power will then shut down automatically.

Warning

Do not use the body fat reading feature of this product if you have a pacemaker or other electronic implanted medical device.

Take a Reading

After programming your personal data you are ready to take a reading.

Select Personal Key and Step On

Using your toe, tap your pre-programmed personal key.

The unit will beep again and the display will show "0.0".





Now step onto the platform.

Get Your Readings

Your weight will be shown first. Continue to stand on the platform.

Your Body Fat Percentage (%FAT) and Body Fat Range will appear on the display screen.

Your total Body Water Percentage (%TBW) will appear on the display screen.

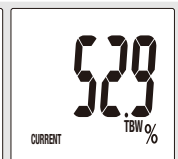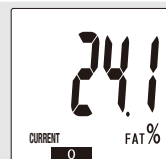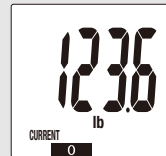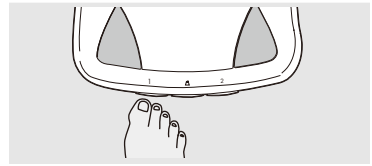Note: Do not step off until Body Fat Percentage is shown.

To obtain previous readings, press Set button while the readings are being displayed.

Press Set button again to go back to the current result display.

Understanding Your Readings

Healthy Range Indicator

Following your body fat percentage reading, a black bar will flash along the bottom of the display, identifying where you fall within the Body Fat Ranges for your age and gender.

Underfat Healthy Overfat Obese


Note: Range of user age is between 7-99

Guest Mode

Guest mode allows you to programme the unit for a one-time use without resetting a personal data number.


Press the Up button to activate the unit.

Then follow the steps: set age, select female or male and specify height.

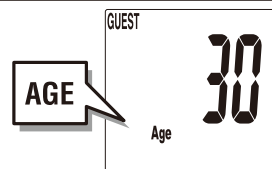Select Weight-Only Mode

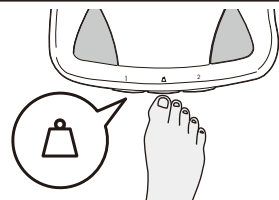Tap the Weight-Only key to select Weight-only mode.

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TANITA Health Equipment H.K. Ltd.
TANITA Europe B.V.
TANITA (Shanghai) Trading Co., Ltd.
TANITA Corporation of America, Inc.
U.S.A.